



COMPRESSION
MEDICAL
STOCKINGS

GENERAL CATALOGUE
2022



PREVENTIVE
*Soft, Medium and
Strong massage
for daily healthy legs*



THERAPEUTIC
*Medical and
anti-embolism stockings
for therapy*



ORTHO
Orthopedic braces

 MADE IN ITALY

Kàmila Compression

Since their introduction, medical compression stockings constituted the treatment having the largest impact on venous stasis, becoming soon the most preferred conservative approach.

They not only apply an external counterpressure that helps a better bloodflow to the heart, resolving thus stasis, but also bring benefits to the superficial and deep venous system, to the blood volume, to the tissues, to the microcirculation and to possible thrombosis.

The compression of elastic stockings is expressed in millimetres of Mercury (mmHg) and it is 100% of the nominal value at the ankle, 70% below the knee and 40% at the thigh.

Compression stockings are divided into two main categories: **PREVENTIVE** and **THERAPEUTIC**.

Within therapeutic stockings, the pressure generated by the compression has to remain constant for 6 months from the first day of use, as long as the instructions for use are followed. The therapeutic stocking has to be well tolerated by the skin of the user, therefore it should not contain harmful yarns or substances, whereas elasticity should not decrease during the first 6 months of use.



RAL CERTIFICATE:
Guaranteed efficacy!

The RAL certificate was introduced in Germany in 1955 by the German Society of Phlebology, establishing the parameters of quality, durability and compression to be met by medical stockings. A German independent body for certification carries out controls and grants RAL certificate only to products in line with their strict parameters. Kamila therapeutic stockings carry RAL certificate and are therefore guaranteed.





COMPANY PROFILE

CA-MI is an Italian family-run factory founded in the early 80's and committed in the production of electromedical equipments, such as SURGICAL SUCTION UNITS for the aspiration of human and animal fluids, used by home-care, hospital, emergency and aesthetic professionals, AEROSOL DELIVERY SYSTEMS used in aerosoltherapy for the treatment of respiratory pathologies, and INHALATORS for hot-humid therapies in the upper airways management. Other medical devices such as diagnostics, body temperature management, TENS devices and antidecubitus products complete the vast range.

The Kàmila brand is created in 2010, as a new line of medical compression stockings, all of them made in Italy. Developed by a management with a long experience in the field, it's a fresh and young product line, with high quality standards. Fine raw materials, state-of-the-art machinery, strict quality controls are the principles of Kàmila. It was on such basis that two lines have been designed: the SUPPORT focusing on prevention, and the MEDICAL dedicated to therapeutic products.

The SUPPORT line pays great attention to fashion by constantly designing new, elegant and refined models. Prevention to be happily lived as an aesthetic pleasure.

The MEDICAL line, carrying the RAL Quality Certificate, offers different medical devices for the management of vascular pathologies, helping the restoration and maintenance of healthy legs, while providing great discretion and comfort at the same time.

CA-MI and KAMILA are always at customer's disposal, counting on in-house R&D Department, and Customer Service.

We welcome our customers and visitors at the major international medical exhibitions, such as MEDICA (Dusseldorf), ARAB HEALTH (Dubai), FIME (USA) and others.

CA-MI manufactures in Italy complying with the UNI EN ISO 9001:2008 and EN ISO 13485:2012/AC:2012 normatives and carry relevant CE 0123 number certified by German Notified Body TÜV SÜD Product Service GmbH.

Kàmila Compression

The claimed compression value is nominal and gradual along the leg. Actually, the indicated compression is applied at the ankles, then it decreases along calves and is lower on thighs. Graduated compression may vary to suit the users' needs. Whoever simply seeks a distressing effect will choose a light compression, whereas firmer compression is required to provide more effective treatment or in presence of pathologies requiring a therapy.

THIGH EDGE
40%
compression

CALF
70%
compression

ANKLE
100%
compression



COMPRESSION LEVELS PREVENTIVE and THERAPEUTIC

PREVENTIVE

Deniers (weight of yarn)	Compression mmHg	Action	Benefits
40	9 - 11	Light Massage	Provides relief from tired legs.
70 70 plus	12 - 15 13 - 16	Medium Massage	Provides relief from tired legs and prevents swelling.
140	18 - 20	Strong Massage	Provides relief from tired legs and prevents swelling. De-stressing.

THERAPEUTIC

	Compression	Indications
KKL I	18 - 21 mmHg	<ul style="list-style-type: none"> ▪ Heaviness and fatigue of the legs ▪ Mild varicose veins during pregnancy ▪ Mild varicose veins without significant oedema ▪ Prophylaxis of thrombosis and embolism for immobile patients ▪ ECS prophylaxis (Economy Class Syndrome during travel)
KKL II	23 - 32 mmHg	<ul style="list-style-type: none"> ▪ More severe varicose veins during pregnancy ▪ Varicose veins with mild oedema ▪ Post superficial and deep venous thrombosis ▪ Post sclerotherapy, venous surgery or thermal ablation ▪ Chronic Venous Insufficiency (CVI) ▪ Management of healed venous ulcer ▪ Prevention of DVT and oedema ▪ Post-thrombotic syndrome (PTS) ▪ Prevention of PTS ▪ Post-traumatic oedema ▪ Lymphoedema
CLASS 2	30 - 40 mmHg	<ul style="list-style-type: none"> ▪ Pronounced varicosis with severe oedema ▪ Severe chronic venous insufficiency ▪ Treatment of active ulcers and recurrent ulcers ▪ Severe post-traumatic oedema ▪ Lymphoedema

DESIGNS & COLOURS

DESIGNS



Knee-high or A-D with closed toe
Knee-high or A-D with open toe



Stay-up or A-G with closed toe
Stay-up or A-G with open toe



Pantyhose or A-T with closed toe
Pantyhose or A-T with open toe



Maternity pantyhose or ATG with closed toe
Maternity pantyhose or ATG with open toe



Thigh (One-Leg) or AGT with closed toe
Thigh (One-Leg) or AGT with open toe

COLOURS



HONEY



WHITE



PLAYA



BEIGE



MINK



BLACK



BLUE



BLACK

INDEX

KNEE-HIGHS PREVENTIVE	<p>LADIES <small>Knee-highs</small></p> <ul style="list-style-type: none"> BAS NEW 70 BAS EASY 70 BAS NEW 140 BAS EASY 140 	<p>UNISEX <small>Knee-highs</small></p> <ul style="list-style-type: none"> FANTASY 15 LONDON 15 COTTON 15 TRAVELLING 15 COTTON 18 	<i>page 8 9</i>
KNEE-HIGHS THERAPEUTIC	<ul style="list-style-type: none"> 7001 A-D 2001 A-D 5001 A-D 5ME1 A-D 5FS1 A-D 	<ul style="list-style-type: none"> 7002 A-D 2002 A-D 3002 A-D 1001 A-D KA A-D 	<i>page 10 11</i>
STAY-UP STOCKINGS PREVENTIVE	<ul style="list-style-type: none"> BORDO COSCIA 70 UP 70 IRIS 70 UP EASY 70 	<ul style="list-style-type: none"> BORDO COSCIA 140 UP 140 IRIS 140 UP EASY 140 	<i>page 12 13</i>
STAY-UP STOCKINGS THERAPEUTIC	<ul style="list-style-type: none"> 7001 A-G 2001 A-G 5001 A-G 	<ul style="list-style-type: none"> 7002 A-G 2002 A-G 3002 A-G <ul style="list-style-type: none"> ■ KA A-G 	<i>page 14 15</i>
PANTYHOSE PREVENTIVE	<ul style="list-style-type: none"> ■ NUDE 40 ■ PRETTY 70 ■ NUDE 70 ■ MICRO 70 ■ GLAMOUR 70 	<ul style="list-style-type: none"> ■ CHERIE 70 ■ FASHION 70 ■ COMFORT 70 ■ COMFORT EASY 70 ■ OVER 70 	<i>page 16 17</i>
PANTYHOSE PREVENTIVE	<ul style="list-style-type: none"> ■ MATERNITY 140 ■ MICRO LADY 140 ■ MABEL 140 	<ul style="list-style-type: none"> ■ COMFORT 140 ■ COMFORT EASY 140 ■ OVER 140 	<i>page 18 19</i>
PANTYHOSE THERAPEUTIC	<ul style="list-style-type: none"> ■ 7001 A-T ■ 2001 A-T ■ 7002 A-T ■ 2002 A-T 	<ul style="list-style-type: none"> ■ 2001 ATG PREGNANCY ■ 2002 ATG PREGNANCY 	<i>page 20 21</i>
THIGH - ONE LEG THERAPEUTIC	<ul style="list-style-type: none"> ■ 7001 AGT ■ 2001 AGT ■ 5001 AGT 	<ul style="list-style-type: none"> ■ 7002 AGT ■ 2002 AGT ■ 3002 AGT <ul style="list-style-type: none"> ■ KA AGT 	<i>page 22 23</i>
THERAPEUTIC KIT	<ul style="list-style-type: none"> ■ ULCER SAN POST SURGERY 		<i>page 24 25</i>
ACCESSORIES THERAPEUTIC	<ul style="list-style-type: none"> ■ LINER FOOT SLIP 	<ul style="list-style-type: none"> ■ SINGLE-USE SOCKS 	<i>page 26 27</i>
USEFUL TIPS SIZE TABLES	<ul style="list-style-type: none"> ■ USEFUL TIPS, MEASUREMENTS POINTS, SIZE TABLES 		<i>page 28 29</i>
ORTHOPEDIC BRACES	<ul style="list-style-type: none"> ■ ANKLE, WRIST, KNEE and GIRDLES 		<i>page 30</i>
Kàmila MERCHANDISING	<ul style="list-style-type: none"> ■ EXHIBITORS, DISPLAYS, ACCESSORIES and GIVEAWAYS 		<i>page 31</i>

KNEE-HIGHS

PREVENTION for LADIES and MEN

Compression Level	Den	mmHg	Product	Designed for	Toe	Sizes	Colours	Composition	Model
MEDIUM MESSAGE Relief from tired legs. Prevent swelling.	70	12 - 15	BAS NEW 70	Ladies	Closed	1 2 3 4 5	Black Playa Black Mink	83% Polyamide 17% Elastam	
MEDIUM MESSAGE Relief from tired legs. Prevent swelling.	70	12 - 15	BAS EASY 70	Ladies	Open	2 3 4 5	Playa Black	83% Polyamide 17% Elastam	
STRONG MESSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	18 - 20	BAS NEW 140	Ladies	Closed	1 2 3 4 5	Black Playa Black Mink	75% Polyamide 25% Elastam	
STRONG MESSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	18 - 20	BAS EASY 140	Ladies	Open	2 3 4 5	Playa Black	79% Polyamide 21% Elastam	
MEDIUM MESSAGE Relief from tired legs. Prevent swelling. Soft cotton on the skin.		15 - 18	COTTON 15	Unisex Travel	Closed	S M L XL	Blue Black Grey	55% Polyamide 32% Cotton 13% Elastam	
MEDIUM MESSAGE Relief from tired legs. Prevent swelling.		14 - 16	FANTASY 15	Unisex Travel	Closed	S M L XL	Blue with bluette lines Black with grey dots	92% Polyamide 8% Elastam	
MEDIUM MESSAGE Relief from tired legs. Prevent swelling.		14 - 16	LONDON 15	Unisex Travel	Closed	S M L XL	Blue with light blue and lavender rhombuses Black with grey and burgundy rhombuses	92% Polyamide 8% Elastam	
MEDIUM MESSAGE Relief from tired legs. Prevent swelling. Soft cotton on the skin and Qskin® microfibre.		15 - 18	TRAVELLING 15	Unisex Travel	Closed	S M L XL XXL	Blue Black Grey	57% Polyamide 33% Cotton 10% Elastam	
STRONG MESSAGE Relief from tired legs. Prevent swelling. Destressing. Soft cotton on the skin and Qskin® microfibre to ensure high breathability and comfort.		18 - 20	COTTON 18	Unisex ravel	Closed	S M L XL XXL	Blue Black Grey	55% Polyamide 33% Cotton 12% Elastam	

LADIES KNEE-HIGHS

	1	2	3	4	5
HEIGHT	155 / 160	160 / 170	170 / 175	175 / 180	180 / 185
SHOE SIZE	35 / 36	36 / 37	38 / 39	39 / 40	41 / 42

UNISEX and MEN KNEE HIGHS

	small	medium	large	x large	xx large
HEIGHT	160 / 170	170 / 180	180 / 185	185 / 190	190 / 195
SHOE SIZE	37 / 39	40 / 42	43 / 44	44 / 45	46 / 47

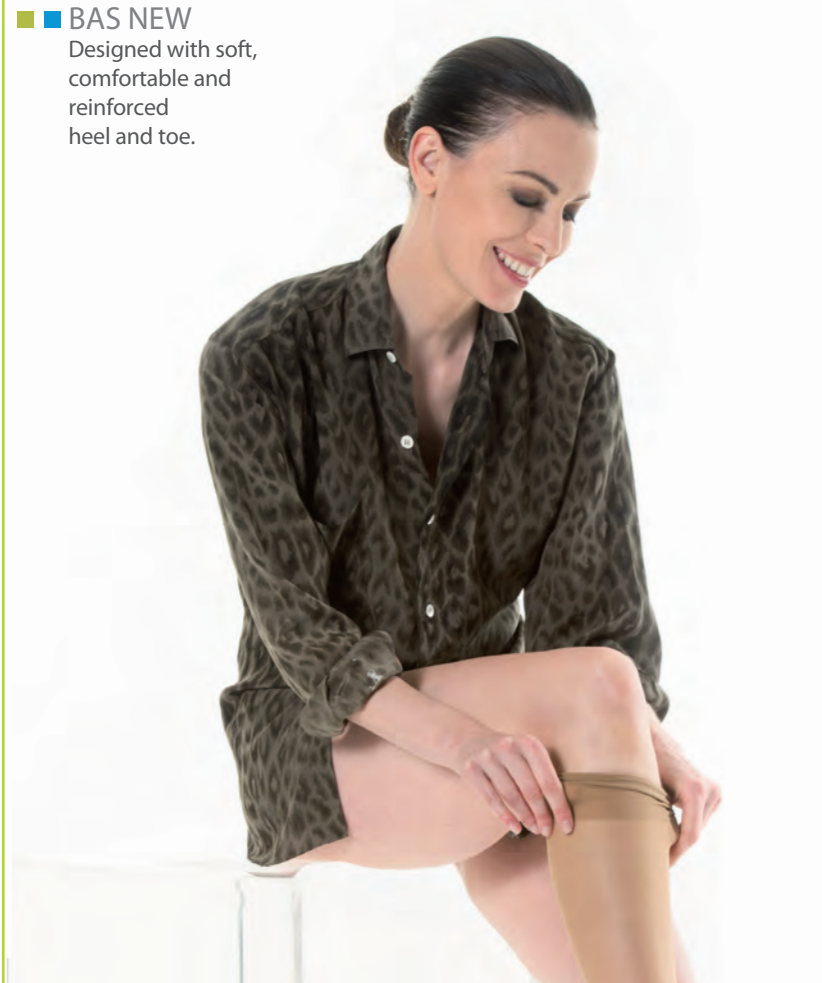
KNEE-HIGHS

PREVENTION for LADIES and MEN



■ ■ **BAS EASY**
 Ideal for shoes with open toe and sandals.
 Comfortable for feet with toe problems.
 Elastic and soft band.

■ ■ **BAS NEW**
 Designed with soft, comfortable and reinforced heel and toe.




■ ■ **COTTON**
 UNISEX
 Indicated for travellers or for those sitting still for many hours.
 Non-constraining elastic band and wide toe for free movements.



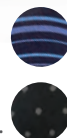■ ■ **LONDON**
 UNISEX
 Characterized by the elegant colourful weave with rhombuses.















■ ■ **FANTASY**
 UNISEX
 With bluette lines or gray polka dots.



KNEE-HIGHS

UNISEX THERAPEUTIC

Compression Level	Product	Toe	Colours	Sizes	Height	Circumf.	Composition	Model	RAL	
KKL I = 18 - 21 mmHg <ul style="list-style-type: none"> ■ Heaviness and fatigue of the legs ■ Mild varicose veins during pregnancy ■ Mild varicose veins without significant oedema ■ Prophylaxis of thrombosis and embolism for immobile patients ■ ECS prophylaxis (Economy Class Syndrome during travel) 	MICRO									
	KT 7001 A-D	Open Closed	Beige Black	S M L XL	■ Regular ■ Long	■ Normal	73% Microfibre 27% Elastam			
	COTY									
	KT 2001 A-D	Open	Beige	S M L XL	■ Regular ■ Long	■ Normal ■ Plus	21% Cotton 53% Polyamide 26% Elastam			
	FUSION									
	FM 5001 A-D	Open	Beige	S M L XL	■ Regular	■ Normal	75% Microfibre 25% Elastam			
FM 5FS1 A-D	Closed	Green	79% Microfibre 21% Elastam							
FM 5ME1 A-D	Closed	Melange (blue)	77% Microfibre 23% Elastam							
KKL II = 23 - 32 mmHg <ul style="list-style-type: none"> ■ More severe varicose veins during pregnancy ■ Varicose veins with mild oedema ■ Post superficial and deep venous thrombosis ■ Post sclerotherapy, venous surgery or thermal ablation ■ Chronic Venous Insufficiency (CVI) ■ Management of healed venous ulcer ■ Prevention of DVT and oedema ■ Post-thrombotic syndrome (PTS) ■ Prevention of PTS ■ Post-traumatic oedema ■ Lymphoedema ■ Lipoedema 	MICRO									
	KT 7002 A-D	Open	Beige Black	S M L XL	■ Regular ■ Long	■ Normal	67% Microfibre 33% Elastam			
	COTY									
	KT 2002 A-D	Open	Beige	S M L XL	■ Regular ■ Long	■ Normal ■ Plus	22% Cotton 49% Polyamide 29% Elastam			
CLASS 2 = 30 - 40 mmHg <ul style="list-style-type: none"> ■ Pronounced varicosis with severe oedema ■ Severe chronic venous insufficiency ■ Treatment of active ulcers and recurrent ulcers ■ Severe post-traumatic oedema ■ Lymphoedema 	COTY ALL USA									
	KU 3002 A-D	Open	Beige	S M L XL	■ Regular ■ Long	■ Normal ■ Plus	31% Cotton 37% Polyamide 32% Elastam			
AFTER-CAST = 15 - 21 mmHg <p>In presence of varicosities, especially after cast removal. Packaged by the piece and not by the pair.</p>	DOGE									
	KT 1001 A-D	Open	Beige	S M L XL XXL	■ Regular	■ Normal	75% Polyamide 25% Elastam			
ANTI-EMBOLISM = 18 - 21mmHg <p>Clinically designed to promote venous return for prevention of thrombosis before, during and after surgery and for medical conditions requiring prolonged bed rest. Before, during and after childbirth.</p>	EMBO-CONTROL									
	KA A-D	Closed with inspection hole	White	S M L XL	■ Regular ■ Long	■ Normal	81% Microfibre 19% Elastam			

KNEE-HIGHS UNISEX THERAPEUTIC



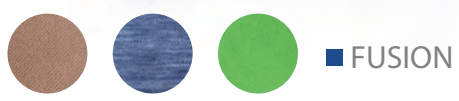
■ MICRO



■ DO-GE



■ EMBO-CONTROL



■ COTY

STAY-UP STOCKINGS

PREVENTION & FASHION

Compression Level	Den	mmHg	Product	Toe	Sizes	Colours	Composition	Model
MEDIUM MASSAGE Relief from tired legs. Prevent swelling.	70	12 - 15	BORDO COSCIA 70	Closed	1 2 3 4 5	Playa Black	88% Polyamide 12% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling.	70	12 - 15	UP 70	Closed	1 2 3 4 5	Playa Mink Black	88% Polyamide 12% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling.	70	12 - 15	IRIS 70	Closed	1 2 3 4 5	Playa Mink Black	76% Polyamide 24% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling.	70	12 - 15	UP EASY 70	Open	2 3 4 5	Playa Black	87% Polyamide 13% Elastam	
STRONG MASSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	18 - 20	BORDO COSCIA 140	Closed	1 2 3 4 5	Playa Black	81% Polyamide 19% Elastam	
STRONG MASSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	18 - 20	UP 140	Closed	1 2 3 4 5	Playa Mink Black	81% Polyamide 19% Elastam	
STRONG MASSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	18 - 20	IRIS 140	Closed	1 2 3 4 5	Playa Mink Black	82% Polyamide 18% Elastam	
STRONG MASSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	18 - 20	UP EASY 140	Open	2 3 4 5	Playa Black	72% Polyamide 28% Elastam	

STAY-UP STOCKINGS

PREVENTION & FASHION



■ ■ UP EASY
Ideal for open toed shoes and sandals and indicated in presence of toe problems. Elastic and soft bands.














■ ■ IRIS
■ ■ UP
Stay-up with 8 cm lace top and double silicone self-adhering band, elegant and handy design with reinforced toe and heel.



■ ■ BORDO COSCIA
Classic stay-ups for garter belts, comfortable and handy, with reinforced toe and heel.

STAY-UP STOCKINGS UNISEX THERAPEUTIC

Compression Level	Product	Toe	Sizes	Height	Circumf.	Colours	Composition	Model	RAL
KKL I = 18 - 21 mmHg <ul style="list-style-type: none"> ■ Heaviness and fatigue of the legs ■ Mild varicose veins during pregnancy ■ Mild varicose veins without significant oedema ■ Prophylaxis of thrombosis and embolism for immobile patients ■ ECS prophylaxis (Economy Class Syndrome during travel) 	MICRO								
	KT 7001 A-G	Open	S M L XL	■ Regular ■ Long	■ Normal	Beige	73% Microfibre 27% Elastam		
	COTY								
	KT 2001 A-G	Open	S M L XL	■ Regular ■ Long	■ Normal ■ Plus	Beige	21% Cotton 53% Polyamide 26% Elastam		
	FUSION								
	FM 5001 A-G	Open	S M L XL	■ Regular	■ Normal	Beige	75% Microfibre 25% Elastam		
KKL II = 23 - 32 mmHg <ul style="list-style-type: none"> ■ More severe varicose veins during pregnancy ■ Varicose veins with mild oedema ■ Post superficial and deep venous thrombosis ■ Post sclerotherapy, venous surgery or thermal ablation ■ Chronic Venous Insufficiency (CVI) ■ Management of healed venous ulcer ■ Prevention of DVT and oedema ■ Post-thrombotic syndrome (PTS) ■ Prevention of PTS ■ Post-traumatic oedema ■ Lymphoedema ■ Lipoedema 	MICRO								
	KT 7002 A-G	Open	S M L XL	■ Regular ■ Long	■ Normal	Beige	67% Microfibre 33% Elastam		
	COTY								
	KT 2002 A-G	Open	S M L XL	■ Regular ■ Long	■ Normal ■ Plus	Beige	22% Cotton 49% Polyamide 29% Elastam		
CLASS 2 = 30 - 40 mmHg <ul style="list-style-type: none"> ■ Pronounced varicosis with severe oedema ■ Severe chronic venous insufficiency ■ Treatment of active ulcers and recurrent ulcers ■ Severe post-traumatic oedema ■ Lymphoedema 	COTY								
	KU 3002 A-G	Open	S M L XL	■ Regular ■ Long	■ Normal ■ Plus	Beige	31% Cotton 37% Polyamide 32% Elastam		
ANTI-EMBOLISM = 18 - 21 mmHg <p>Clinically designed to promote venous return for prevention of thrombosis before, during and after surgery and for medical conditions requiring prolonged bed rest. Before, during and after childbirth.</p>	EMBO-CONTROL								
	KA A-G	Closed with inspection hole	S M L XL	■ Regular	■ Normal ■ Plus	White	81% Microfibre 19% Elastam		

STAY-UP STOCKINGS
UNISEX THERAPEUTIC



■ COTY



■ EMBO-CONTROL



■ MICRO



Silicone "bubbled" bands for better adherence to thigh and sweat drainage.

PANTYHOSE

PREVENTION & FASHION

Compression Level	Den	mmHg	Product	Toe	Sizes	Colours	Composition	Model
LIGHT MASSAGE Relief from tired legs.	40	9 - 11	NUDE 40	Closed	1 2 3 4	Honey Mink Black	80% Polyamide 20% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling.	70	12 - 15	PRETTY 70	Closed	1 2 3 4	Honey Mink Blue Black	79% Polyamide 21% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling.	70	12 - 15	NUDE 70	Closed	1 2 3 4	Honey Mink Black	76% Polyamide 24% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling.	70	12 - 15	MICRO 70	Closed	1 2 3 4 5	Black	76% Polyamide 24% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling.	70	13 - 16	GLAMOUR 70	Closed	1 2 3 4 5	Black	85% Polyamide 15% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling. De-stressing.	70	12 - 15	FASHION 70	Closed	1 2 3 4 5	Mink Playa Black	83% Polyamide 17% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling. De-stressing.	70	12 - 15	COMFORT 70	Closed	1 2 3 4 5*	Playa Mink Black	82% Polyamide 18% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling. De-stressing.	70	12 - 15	COMFORT EASY 70	Open	2 3 4 5*	Playa Black	79% Polyamide 21% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling. De-stressing.	70	12 - 15	OVER 70	Closed	2* 3* 4* 5*	Playa Mink Black	83% Polyamide 17% Elastam	
MEDIUM MASSAGE Relief from tired legs. Prevent swelling.	70	13 - 16	CHERIE 70	Closed	1 2 3 4 5*	Honey Mink Black	87% Polyamide 13% Elastam	

* Plus size

PANTYHOSE

PREVENTION & FASHION

■ FASHION

Up-to-date and comfortable version of the classic pantyhose: soft waistband combined with form-fitting lace panty. Reinforced heel and toe.

■ COMFORT

Classic, long durability, with handy waistband stable for any body shape. Form-fitting lace panty and reinforced heel and toe.

■ COMFORT EASY

Ideal for open toed shoes and sandals and indicated in presence of toe problems. Elastic and soft bands.

■ OVER

Lace panty with special gusset for greater comfort. Ideai to minimise large waistline and hips.

■ PRETTY

Young and refined pantyhose, high leg, embroidered lace panty. Ultra sheer leg and sheer heel.



MICROFIBRE

■ MICRO




Pantyhose in soft microfibre, great wearability. With control top and embroidered high leg cut.

■ NUDE

Sheer to waist and delicate pantyhose, sheer heel and malicious embroidery under the waistband.

PANTYHOSE

PREVENTION & FASHION

Compression Level	Den	mmHg	Product	Toe	Sizes	Colours	Composition	Model
STRONG MESSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	16 - 18	MICRO LADY 140	Closed	1 2 3 4* 5*	Black	79% Polyamide 21% Elastam	
STRONG MESSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	18 - 20	MABEL 140	Closed	1 2 3 4 5*	Honey Playa Mink Black	79% Polyamide 21% Elastam	
STRONG MESSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	18 - 20	COMFORT 140	Closed	1 2 3 4 5*	Playa Mink Black	77% Polyamide 23% Elastam	
STRONG MESSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	18 - 20	COMFORT EASY 140	Open	2 3 4 5*	Playa Black	74% Polyamide 26% Elastam	
STRONG MESSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	18 - 20	OVER 140	Closed	2* 3* 4* 5*	Playa Mink Black	77% Polyamide 23% Elastam	
STRONG MESSAGE Relief from tired legs. Prevent swelling. De-stressing.	140	18 - 20	MATERNITY 140	Closed	2 3 4 5	Playa Black	77% Polyamide 23% Elastam	

* Plus size

PANTYHOSE

PREVENTION & FASHION



■ **CHERIE**
 Perfect blend between refined elegance and relaxing action in innovative tights. Silky knit, sheer heel, embroidery at the waistband.



■ **MABEL**
 Efficient and elegant, with accurate compression. Control top with charming embroidery, reinforced toe for a long durability.



■ **COMFORT EASY**
 Ideal for open toed shoes and sandals and indicated in presence of toe problems. Elastic and soft bands.



■ **OVER**
 Lace panty with special gusset for greater comfort. Ideai to minimise large waistline and hips.



MICROFIBRE

■ **MICRO LADY**
 Tights in soft microfibre, great wearability. With control top and embroidered high leg cut.

■ **COMFORT**
 Classic, long durability, with handy waistband stable for any body shape. Form-fitting lace panty and reinforced heel and toe.

PANTYHOSE

UNISEX & PREGNANCY THERAPEUTIC

Compression Level	Product	Toe	Sizes	Height	Circumf.	Colours	Composition	Model	RAL
KKL I = 18 - 21 mmHg <ul style="list-style-type: none"> ■ Heaviness and fatigue of the legs ■ Mild varicose veins during pregnancy ■ Mild varicose veins without significant oedema ■ Prophylaxis of thrombosis and embolism for immobile patients ■ ECS prophylaxis (Economy Class Syndrome during travel) 	MICRO								
	KT 7001 A-T	Open	S M L XL	■ Regular ■ Long	■ Normal	Beige	73% Microfibre 27% Elastam		
	COTY								
	KT 2001 A-T	Open	S M L XL	■ Regular ■ Long	■ Normal	Beige	21% Cotton 53% Polyamide 26% Elastam		
	COTY								
	KT 2001 ATG PREGNANCY	Closed	S M L XL	■ Regular ■ Long	■ Normal	Beige	21% Cotton 53% Polyamide 26% Elastam		
KKL II = 23 - 32 mmHg <ul style="list-style-type: none"> ■ More severe varicose veins during pregnancy ■ Varicose veins with mild oedema ■ Post superficial and deep venous thrombosis ■ Post sclerotherapy, venous surgery or thermal ablation ■ Chronic Venous Insufficiency (CVI) ■ Management of healed venous ulcer ■ Prevention of DVT and oedema ■ Post-thrombotic syndrome (PTS) ■ Prevention of PTS ■ Post-traumatic oedema ■ Lymphoedema ■ Lipooedema 	MICRO								
	KT 7002 A-T	Open	S M L XL	■ Regular ■ Long	■ Normal	Beige	67% Microfibre 33% Elastam		
	COTY								
	KT 2002 A-T	Open	S M L XL	■ Regular ■ Long	■ Normal	Beige	22% Cotton 49% Polyamide 29% Elastam		
	COTY								
	KT 2002 ATG PREGNANCY	Open	S M L XL	■ Regular ■ Long	■ Normal	Beige	22% Cotton 49% Polyamide 29% Elastam		

PANTYHOSE

UNISEX & PREGNANCY THERAPEUTIC



■ PREGNANCY














■ MICRO
■ COTY

THIGH - ONE LEG

THERAPEUTIC

Left, Right or Ambidextrous

Compression Level	Product	Toe	Sizes	Height	Circumf.	Colours	Composition	Unit	RAL	
KKL I = 18 - 21 mmHg <ul style="list-style-type: none"> ■ Heaviness and fatigue of the legs ■ Mild varicose veins during pregnancy ■ Mild varicose veins without significant oedema ■ Prophylaxis of thrombosis and embolism for immobile patients ■ ECS prophylaxis (Economy Class Syndrome during travel) 	MICRO									
	KT 7001 AGT Left, Right or Ambidextrous	Open Closed	S M L XL	■ Regular ■ Long	■ Normal	Beige Black	73% Microfibre 27% Elastam			
	COTY									
	KT 2001 AGT Left, Right or Ambidextrous	Open	S M L XL	■ Regular ■ Long	■ Normal ■ Plus	Beige	21% Cotton 53% Polyamide 26% Elastam			
	FUSION									
	FM 5001 AGT Left, Right or Ambidextrous	Open	S M L XL	■ Regular	■ Normal	Beige	75% Microfibre 25% Elastam			
	KKL II = 23 - 32 mmHg <ul style="list-style-type: none"> ■ More severe varicose veins during pregnancy ■ Varicose veins with mild oedema ■ Post superficial and deep venous thrombosis ■ Post sclerotherapy, venous surgery or thermal ablation ■ Chronic Venous Insufficiency (CVI) ■ Management of healed venous ulcer ■ Prevention of DVT and oedema ■ Post-thrombotic syndrome (PTS) ■ Prevention of PTS ■ Post-traumatic oedema ■ Lymphoedema ■ Lipoedema 	MICRO								
		KT 7002 AGT Left, Right or Ambidextrous	Open Closed	S M L XL	■ Regular ■ Long	■ Normal	Beige Black	67% Microfibre 33% Elastam		
	COTY									
	KT 2002 AGT Left, Right or Ambidextrous	Open	S M L XL	■ Regular ■ Long	■ Normal ■ Plus	Beige	22% Cotton 49% Polyamide 29% Elastam			
CLASS 2 = 30 - 40 mmHg <ul style="list-style-type: none"> ■ Pronounced varicosity with severe oedema ■ Severe chronic venous insufficiency ■ Treatment of active ulcers and recurrent ulcers ■ Severe post-traumatic oedema ■ Lymphoedema 	COTY ALL USA									
	KU 3002 AGT Left or Right	Open	S M L XL	■ Regular ■ Long	■ Normal ■ Plus	Beige	31% Cotton 37% Polyamide 32% Elastam			
ANTI-EMBOLISM = 18 - 21 mmHg <p>“Embo-Control” feature anatomic design and great wearability. All designs have inspection hole useful in the prevention of mycosis, and anatomic shape anti-irritation heel pocket. Clinically designed to promote venous return for prevention of thrombosis before, during and after surgery and for medical conditions requiring prolonged bed rest. Before, during and after childbirth.</p>	EMBO-CONTROL									
	KA AGT Left, Right or Ambidextrous	Closed with inspection hole	S M L XL	■ Regular ■ Long	■ Normal ■ Plus	White	81% Microfibre 19% Elastam			

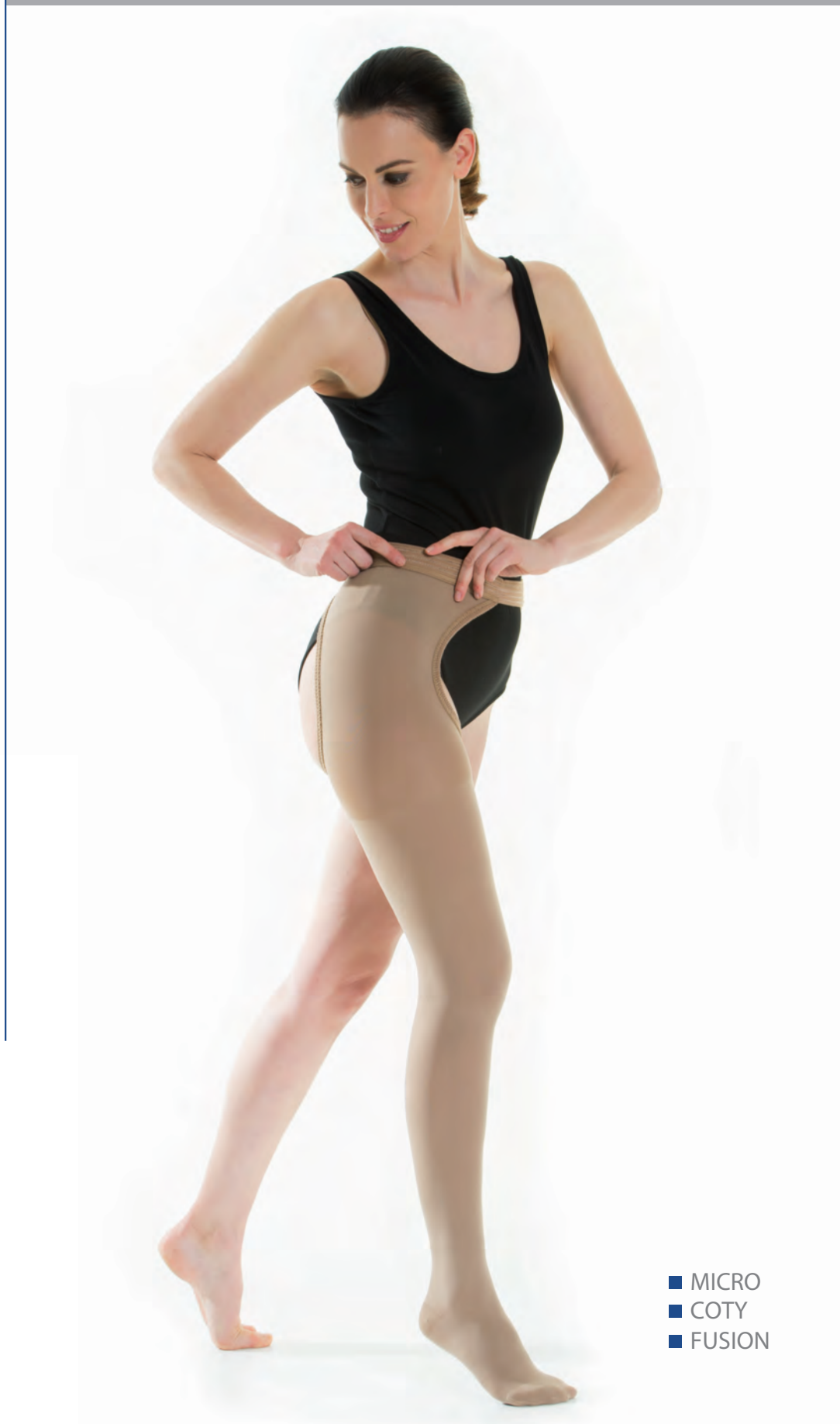
THIGH - ONE LEG
THERAPEUTIC
Left, Right or Ambidextrous



■ **EMBO-CONTROL**
Thigh (One-Leg) with
adjusting belt,
ambidextrous.



■ **EMBO-CONTROL**
Inspection hole



■ **MICRO**
■ **COTY**
■ **FUSION**

KNEE-HIGHS | THIGH, ONE-LEG THERAPEUTIC KITS

ULCER SAN - Available with two different compression levels, the ULCER SAN kits combine a therapeutic cotton knee-high with an anti-embolism knee-high, constituting an advanced therapy as an efficient alternative to multiple bendages, indicated for the treatment of cutaneous ulcers with or without venous insufficiency.

Product	Total Compression	Toe	Sizes	Height	Circumf.	Colours	Composition
ULCER-SAN 2001							
US 2001 A-D 1 x COTY Knee-High Compression 18 - 21 mmHg	36 - 39 mmHg	Open	S M L XL	▪ Regular	▪ Normal	Beige	21% Cotone 53% Polyamide 26% Elastam
KA A-D 1 x Embo Control Knee-High Compression 18 - 21 mmHg		Closed with inspection hole	S M L XL	▪ Regular	▪ Normal	White	81% Microfibre 19% Elastam
KM 026 2 x A-D Liners		Closed	S/M L/XL	-	-	White	100% Polyamide
KM 017 1 x Foot Slip		-	-	-	-	White	100% Polyamide
ULCER-SAN 2002							
US 2002 A-D 1 x COTY Knee-High Compression 23 - 32 mmHg	41 - 45 mmHg	Open	S M L XL	▪ Regular	▪ Normal	Beige	22% Cotone 49% Polyamide 29% Elastam
KA A-D 1 x Embo Control Knee-High Compression 18 - 21 mmHg		Closed with inspection hole	S M L XL	▪ Regular	▪ Normal	White	81% Microfibre 19% Elastam
KM 026 2 x A-D Liners		Closed	S/M L/XL	-	-	White	100% Polyamide
KM 017 1 x Foot Slip		-	-	-	-	White	100% Polyamide

POST SURGERY - The POST-SURGERY kit includes a special Thigh (One-Leg) made with cotton in reduced thickness providing a great fit, even wearing regular shoes. This particular Thigh (One-Leg) has a compression of 23-32 mmHg at the ankle, gradually decreasing to 18 mmHg at the thigh, allowing the compliance of the prescribed therapy after surgery.

Product	Compression	Toe	Sizes	Height	Circumf.	Colours	Composition
POST-SURGERY							
PS 2002 1 x 2002 AGT Thigh (One-Leg) Left or Right	▪ 23 -32 mmHg at the ankle ▪ 18 mmHg at the thigh	Open	S M L XL	▪ Regular ▪ Long	▪ Plus	Beige	24% Cotone 43% Polyamide 33% Elastam
KM 028 2 x A-G Liner		Closed	S/M L/XL	-	-	White	100% Polyamide
KM 017 1 x Foot Slip		-	-	-	-	White	100% Polyamide

KNEE-HIGHS | THIGH, ONE-LEG
THERAPEUTIC KITS



■ ULCER-SAN




■ POST-SURGERY

ACCESSORIES THERAPEUTIC

Description	Product	Toe	Sizes	Colours	Composition	Unit
LINER						
To be worn under the stocking to support and protect dressings and the skin, to aid wearing stockings.	KM 026 	Closed	S / M L / XL	White	100% Polyamide	By pair
	KM 028 	Closed	S / M L / XL	White	100% Polyamide	By pair
FOOT SLIP						
To aid wearing stockings.	KM 017	Closed	One size fits all	White	100% Polyamide	By piece
SINGLE-USE SOCKS						
Sanitary Disposable Try on Socks. Supplied in convenient dispenser box of 50's which is easily stored in changing rooms, seating areas or at till points for staff to issue. Box is perforated for easy, neat opening. Perfect for shoe shops, department stores, boutiques, market stalls, sports shops or anywhere selling shoes, boots and other footwear.	KM 055	Closed	One size fits all	White with assorted coloured circles	100% Polyamide	By pair in box



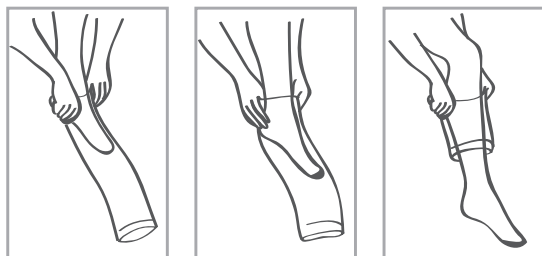HOW TO WEAR CLOSED TOE STOCKINGS



- Unroll the stocking down to the heel and wear it evenly over your foot.
- Pull up the stocking evenly without tugging. Spread the fabric over the leg using the flat of your hand.
- Continue to pull it up until it perfectly fits all over your leg.

HOW TO WEAR OPEN TOE STOCKINGS



- Place your foot into the foot slip.
- Put the stocking on.
- Remove the foot slip by pulling its tip.

UNIVERSAL APPLICATOR FOR STOCKINGS

EZY-AS

Description	Product	Sizes	Colours	Circumferences A or B
APPLICATOR EZY-AS is a universal applicator suitable for stockings with open toe as well as with closed toe, too. EZY-AS comes in three colour coded sizes to suit, each one with smooth lines, no corners to snag garments and can be used by patients on their own or by carers of patients who are unable to apply their own garments. It is engineered for durability, it is lightweight, guaranteed and it can be cleaned with soap and water or hypochlorites if being used in multiple patient use areas to prevent HAI. EZY-AS design means it will not rust, chip, peel, leak, snag or damage garments. Ready to use out of the box.	SAS 08017	S	Red ■	28 – 35 cm
	SAM 06264	M	Yellow ■	32 – 40 cm
	SAL 06264	L	Blue ■	38 – 47 cm



Tutorial video available upon request

HANDLES EZY-AS offers the option of a smart universal handle system which fits all three sizes and can be easily removed using the integral clips when not required. This enables extendable reach for those with limited movement or require this during their recuperation time from surgery.	UHA 07246	One size fits all	Black ■	
-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------	-----------	-------------------	---------	--



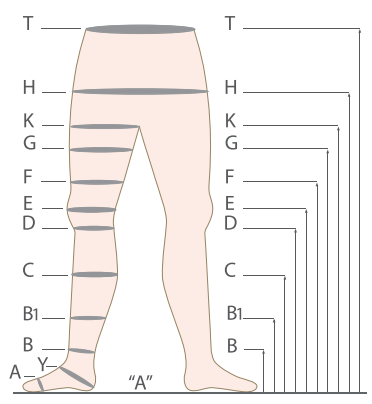
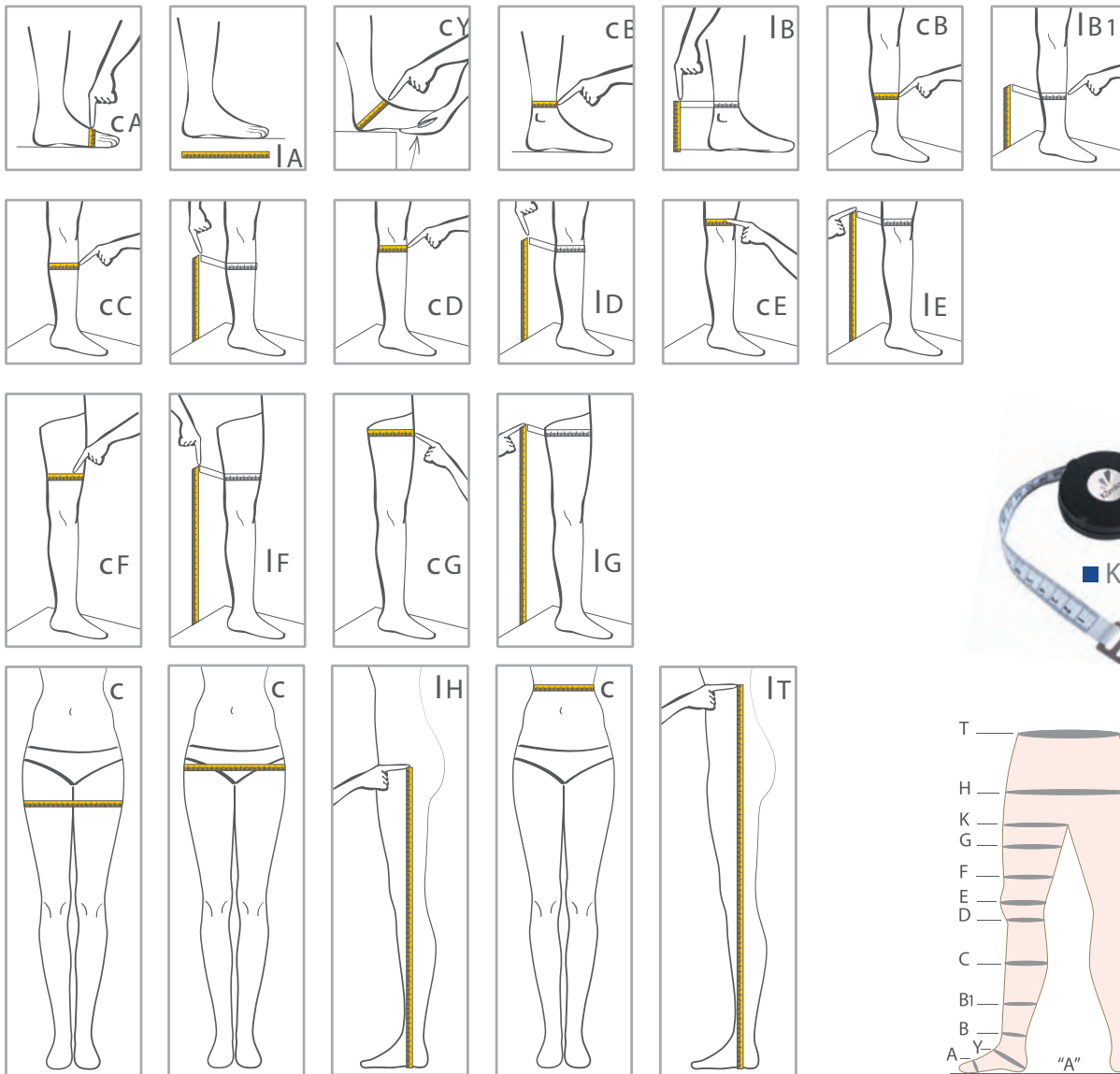
INSTRUCTIONS

Position the applicator on a firm surface at a comfortable working height with the elongated hole facing you. Insert the stocking into the applicator, ensuring the heel of the stocking is centered as you lower the garment into the applicator. Place the applicator on the floor and lower the foot into the stocking. Pull the applicator upwards, up to the knee as indicated in the pictures.



HOW TO WEAR USEFUL TIPS

- Fingernails, rings, bracelets, etc. may accidentally damage the Kàmila Stockings. Remove them or use gloves to protect the stockings.
- For those who have limited hand strength the use of extra grip gloves can help tremendously.
- Stockings should be washed everyday by hand with neutral detergents, temperature should not exceed 30°C (86°F).
- To dry the stockings lay them on a flat surface (do not hang), such as a spongy towel and keep away from heat sources or direct light.
- Check your shoes carefully. Eliminate any rough or sharp surfaces that may damage the stockings.
- Any loose threads inside the heel are not defects and shall not be cut off.
- Stockings should be put on upon rising in the morning and taken off before going to bed.
- Should you have to wear the stockings for a long period, an additional pair is recommended for daily alternating. This will provide optimal results and therapeutic efficacy.
- The Kàmila Stockings should be stored in a place with temperature between 5°C and 25°C (41°F and 77°F).



SIZE TABLES

THERAPEUTIC, KITS and EMBO-CONTROL

MICRO | FUSION | ULCER SAN

		small	medium	large	xlarge
cm - circumference	cT	67-70	71-80	81-90	91-100
	cH	85-90	91-100	101-110	111-120
	cG	46-50	52-57	58-63	64-70
	cD	29-32	34-37	38-42	43-47
	cC	31-34	36-39	40-44	45-49
	cB	20-22	23-25	26-28	29-31

COTY

		small	small plus	medium	medium plus	large	large plus	xlarge	xlarge plus
cm - circumference	cT	67-70		71-80		81-90		91-100	
	cH	85-90		91-100		101-110		111-120	
	cG	46-50	46-55	52-57	52-62	58-63	58-69	64-70	64-77
	cD	29-32	29-37	34-37	34-42	38-42	38-48	43-47	43-53
	cC	31-34		36-39		40-44		45-49	
	cB	20-22		23-25		26-28		29-31	

DOGE

		small	medium	large	xlarge	xxlarge
cm circumf.	cC	30-33	33-38	38-42	42-45	45-48
	cB	19-21	22-24	25-27	28-29	30-31

POST SURGERY

		small plus	medium plus	large plus	xlarge plus
cm - circonf.	cT	67 - 70	71 - 80	81 - 90	91 - 100
	cG	46 - 55	52 - 62	58 - 69	64 - 77
	cD	29 - 37	34 - 42	38 - 48	43 - 53
	cC	31 - 34	36 - 39	40 - 44	45 - 49
	cB	20 - 22	23 - 25	26 - 28	29 - 31




EMBO CONTROL

HEIGHT		
ARTICLE	TYPE	ARTICLE SIZE height to point to point
A-D KNEE -HIGH	regular	from A to ℓD < 39 cm
A-D KNEE -HIGH	long	from A to ℓD > 39 cm
A-G STOCKING AGT THIGH (one leg)	regular	from A to ℓG < 72 cm
A-G STOCKING AGT THIGH (one leg)	long	from A to ℓG > 72 cm

CIRCUMFERENCE					
	circumf. C (cm)	small	medium	large	xlarge
ALL MODELS	c B	18-21	22-24	25-27	28-29
ALL MODELS	c C	30-33	33-38	38-42	42-45
A-G STOCKING ambidextrous AGT THIGH (one leg) ambidextrous	c G	< 48	< 55	< 62	< 74
A-G STOCKING plus AGT THIGH plus - dx or sn	c G	< 55	< 62	< 69	< 76

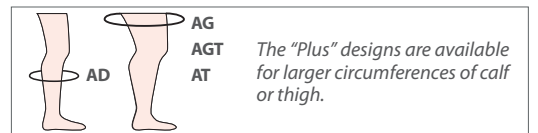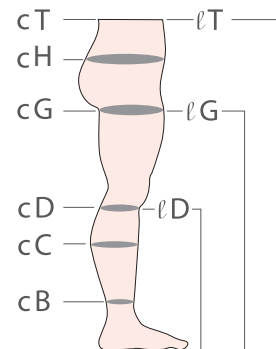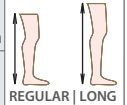
Reference points for measurements

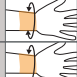

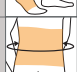


C = circumference ℓ = height	
cT	hips
cH	pelvis (largest point)
cG	5 cm below groin
cD	3 cm below knee pit (Popliteal fossa)
cC	at mid calf
cB	ankle over malleolus (tightest point)
ℓT	from A heel baseline to cT
ℓG	from A heel baseline to cG REGULAR: Up to 72 cm / LONG above 72 cm
ℓD	from A heel baseline to cD REGULAR: Up to 39 cm / LONG above 39cm



ORTHOPEDIC BRACES

ANKLE, WRIST, KNEE and GIRDLES

Characteristics	Product	Sizes	Colours	Composition	
ANKLE SUPPORTS					
The anatomical shaped knee supports in knitted fabric are indicated to support joints and to recover during the post-trauma phase, providing great comfort thanks to its bi-expandable fabric and flat seams. Both models have the same characteristics.	8044-O Ortho Ankle	S M L XL	Beige	47% Cotton 36% Polyamide 17% Elastam	
	8044-S Sport Ankle	S M L XL	Blue White		
WRIST BRACES					
The elastic wrist braces protect joints during sport practice. Indicated for prevention of sprains and joint pains, and post-cast. The anatomical shape does not prevent normal movements and the cotton provides great comfort. Height 8 cm.	8022-O Ortho Wrist	S M L	Beige White	69% Cotton 20% Polyamide 11% Elastam	
	8022-S Sport Wrist	S M L XL	White with coloured stripes		
KNEE SUPPORTS					
The Ortho, Sport and Winter knee supports protect joints during sport practice, after trauma and post-cast. The anatomical shape does not prevent normal movements and all the three models provide great comfort. Both Ortho and Sport have the same characteristics, whereas Winter in wool is indicated for rheumatism.	8062-O Ortho Knee	S M L XL	Beige	49% Cotton 20% Viscosa 16% Polyamide 11% Rubber 4% Elastam	
	8062-S Sport Knee	S M L XL	Blue White		
	8062-L Winter Knee	S M L XL	Beige		
GIRDLES					
The abdominal girdles in cotton with velcro closure (single or double) are indicated for the post-operative treatment and prophylaxis, as well as for the prophylaxis post-partum when the pressure exerted from the abdominal wall has to be reduced. Girdles height: 27 cm. Great comfort of cotton on the skin.	9043 Single Velcro	S M L XL XXL XXXL XXXXL	White	73% Cotton 15% Polyamide 12% Rubber	
	9044 Double Velcro				

PRODUCT		small	medium	large	xlarge	xxlarge	xxxlarge	xxxxlarge
Ortho Wrist cm / wrist circumference		16 - 21	19 - 25	23 - 29	//	//	//	//
Sport Wrist cm / wrist circumference		16 - 21	19 - 25	23 - 29	27 - 34	//	//	//
Ankle Support (all) cm / in the tightest point around the ankle		30 - 35	36 - 39	40 - 44	45 - 49	//	//	//
Knee Support (all) cm / knee circumference		19 - 21	22 - 24	25 - 27	28 - 29	//	//	//
Girdles (all) cm / hips circumference		60 - 70	70 - 85	85 - 100	100 - 115	115 - 130	130 - 145	145 - 160

Kàmila MERCHANDISING

EXHIBITORS, DISPLAYS, ACCESSORIES and GIVEAWAYS

EXHIBITORS

- Ground exhibitor in metal code 165
- Ground exhibitor in cardboard code 151
- Counter exhibitor for pantyhose code 160
- Counter exhibitor for knee-highs code 161



WINDOW DISPLAYS

(Italian language, other languages available upon request)

- Double-face window display code 154
- 65 x 98 cm
- Ezy-As window display code 158
- 45 x 65 cm
- Ezy-As "opentoe" window display code 156
- 45 x 65 cm



ACCESSORIES and GIVEAWAYS



- Colour ring code KM0100
- Measuring tape code KM0200
- Lady leg code KM1419
- Man leg code KM014

Pictures, images and characteristics of the merchandising may vary without notice.



PREVENTION
Soft, Medium and Strong massage for daily healthy legs

THERAPEUTIC
Medical and anti-embolism stockings for therapy

ORTHO
Orthopedic braces



Kàmila by CA-MI srl
Via Ugo La Malfa 13 - Frazione Pilastro - 43013 Langhirano (PR) - Italia
Tel. +39 0521 637133 / 631138 - FAX +39 0521 639041
info@kamilamedical.com - www.kamilamedical.com
www.ca-mi.it
CATALOGUE Kàmila - KM023 Rev. 07 - MARCH 2022

